




SPP

- ▶ **Pre-Ischial Shelf and Lateral Adductor build-up**
- ▶ **Contoured Base, High Density Molded Foam**
- ▶ **Twin-Cell, Aqueous Gel Insert Under At-Risk Bony Areas**



SOLUTION[®] GEL

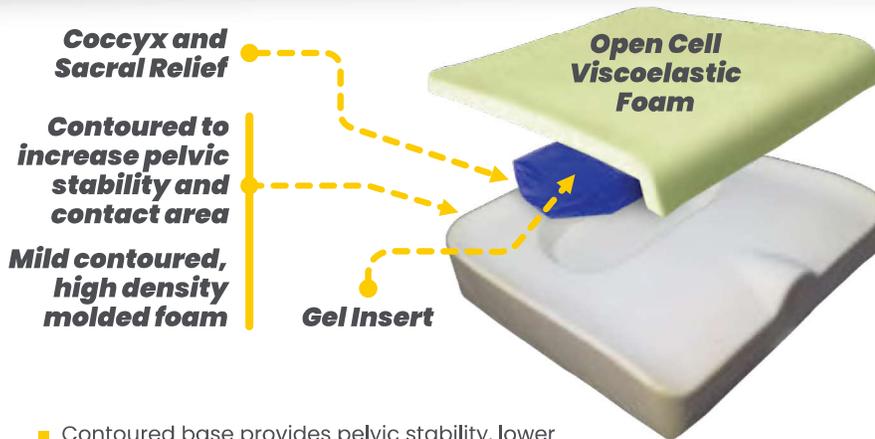

- Contoured base provides pelvic stability, lower extremity support, and increases overall surface contact area.
- Pelvic well provides a pre-ischial shelf to prevent posterior pelvic tilt. Segmented for each IT, an aqueous gel pack provides additional skin protection and reduces shearing forces.
- Top layer of viscoelastic memory foam provides optimal comfort and protection through immersive properties.
- Coolcore is the standard cushion cover material that has proven temperature and moisture control.
- Fire resistance: ISO 8191-1, ISO 8191-2, CAL 117-2013
- Biological evaluation/cytotoxicity ISO 10993-5, 10993-10
- Sub-trochanter support to increase lateral stability
- Maintains alignment of the femurs medially and laterally

WEIGHT CAPACITY

450 lbs.

OVERALL WEIGHT

3½ - 5 lbs.

DIMENSIONS

from 10in (25.40cm) x 10in (25.40cm)
to 22in (55.88cm) x 24in (60.96cm)

SOF¹ available

PROFILE

3 in (7.62cm)

MSRPS

\$375 - \$430

HCPCS

E2607 | E2608

1. Special Order Form
2. CoolCore[®] is standard on all cushions except the Simplicity G.
3. Special order sizes available for all cushions except the Glacial SP. Contact sales for additional assistance
4. Allowances listed are estimations for cushions billed in conjunction with Group 3 power chair bases using the KU modifier and may vary by state