




# SPP

- ▶ **Pre-Ischial Shelf and Lateral Adductor build-up**
- ▶ **Mild Contoured Base, High Density Molded Foam**
- ▶ **Ischial insert consisting of multi-layered viscoelastic foam**



## SPECTRUM FOAM

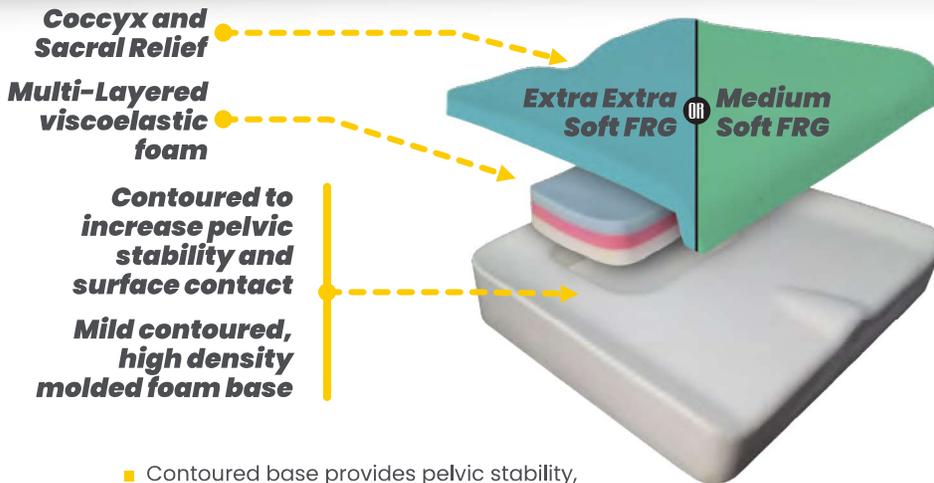

- Contoured base provides pelvic stability, lower extremity support, and increases overall surface contact area
- Pelvic well provides a pre-ischial shelf to prevent posterior pelvic tilt. A triple layer of viscoelastic foam within the well promotes optimal immersion and envelopment
- Top layer of Extra Extra Soft SunMate FRG foam provides optimal protection through immersive properties and improves micro-climate. Medium-Soft SunMate FRG foam used on cushions greater or equal to 22in (55.88cm) wide
- Coolcore is the standard cushion cover material that has proven temperature and moisture control
- Fire resistance: ISO 8191-1, ISO 8191-2, CAL 117-2013
- Biological evaluation/Cytotoxicity ISO 10993-5, 10993-10
- Sub-trochanter support to increase lateral stability
- Maintains alignment of the femurs medially and laterally

### WEIGHT CAPACITY

**650 lbs.**

### OVERALL WEIGHT

**3 lbs.**

### DIMENSIONS

from 14in (35.56cm) x 18in (45.72cm)  
to 24in (60.96cm) x 26in (66.04cm)

SOP<sup>1</sup> available

### PROFILE

**3in (7.62cm)**

### MSRPS

**\$460 - \$590**

### HCPCS

**E2607 | E2608**

1. Special Order Form
2. CoolCore® is standard on all cushions except the Simplicity G.
3. Special order sizes available for all cushions except the Glacial SP. Contact sales for additional assistance
4. Allowances listed are estimations for cushions billed in conjunction with Group 3 power chair bases using the KU modifier and may vary by state