

- ▶ **Extra Soft SunMate FRG Foam Top**
- ▶ **Twin-Cell, Lightweight Fluid Insert Under At-Risk Bony Areas**
- ▶ **Pre-Ischial Shelf**
- ▶ **Mild Contoured Base, High Density Molded Foam**

SPP



SPECTRUM GEL

WEIGHT CAPACITY

550 lbs.

OVERALL WEIGHT

3½ lbs.

DIMENSIONS

from 14in (35.56cm) x 14in (35.56cm)
to 24in (60.96cm) x 24in (60.96cm)

SOF¹ available

PROFILE

3in (7.62cm)

MSRPS

\$460 - \$530

HCPCS


E2607 | E2608

Coccyx and Sacral Relief

Gel Insert

Contoured to increase pelvic stability and surface contact

Mild contoured, high density molded foam base

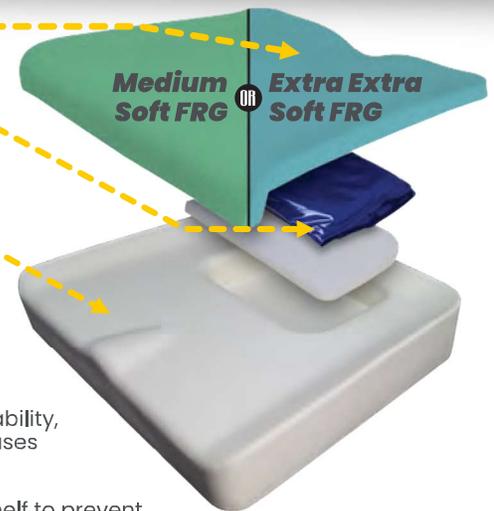

- Contoured base provides pelvic stability, lower extremity support, and increases overall surface contact area
- Pelvic well provides a pre-ischial shelf to prevent posterior pelvic tilt. Segmented for each IT, a lightweight gel pack provides additional skin protection, pelvic stability, and optimal immersion
- Top layer of Extra Extra Soft SunMate FRG foam provides optimal protection through immersive properties and improves micro-climate. Medium Soft SunMate FRG foam used on cushions greater or equal to 22in (55.88cm) wide
- CoolCore is the standard cushion cover material that has proven temperature and moisture control
- Pre-ischial support to help prevent posterior pelvic tilt
- Fire resistance: ISO 8191-1, ISO 8191-2, CAL 117-2013
- Biological evaluation/Cytotoxicity ISO 10993-5, 10993-10
- Sub-trochanter support to increase lateral stability
- Maintains alignment of the femurs medially and laterally

1. Special Order Form
 2. CoolCore® is standard on all cushions except the Simplicity G.
 3. Special order sizes available for all cushions except the Glacial SP. Contact sales for additional assistance
 4. Allowances listed are estimations for cushions billed in conjunction with Group 3 power chair bases using the KU modifier and may vary by state