

- ▶ **Multi-layered foam**
- ▶ **Deep Contour**
- ▶ **Coolcore® cushion cover has pockets on sides**

SPP



TRU-COMFORT® 2

WEIGHT CAPACITY

450 lbs.

OVERALL WEIGHT

4 lbs.

DIMENSIONS

from 14in (35.56cm) x 14in (35.56cm)
to 24in (60.96cm) x 24in (60.96cm)

SOF¹ available

PROFILE

4½ in (11.43cm)

MSRPS

\$400 - \$500


HCPCS

E2607 | E2608

Coccyx and Sacral Relief

Contoured to increase pelvic stability and surface contact

Deeply contoured, high density molded foam base

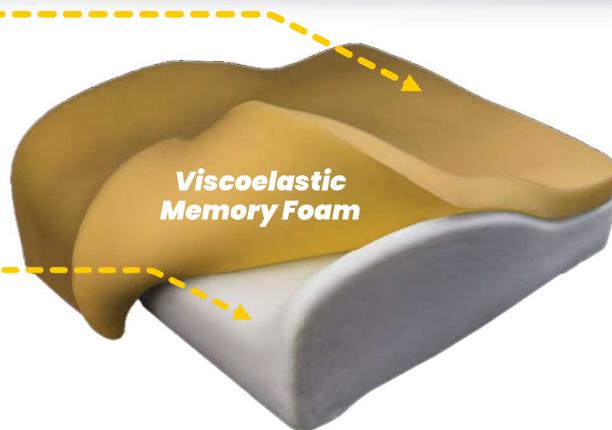

- Deep contoured base provides maximal pelvic stability, lower extremity support, and surface contact area
- Sub-trochanter support to increase lateral stability.
- Maintains alignment of the femurs medially and laterally
- Anti-thrust component provides a pre-ischial shelf to promote pelvic positioning. Abductor and adductor components create troughing to promote lower extremity positioning
- Top layer of viscoelastic memory foam provides optimal comfort and protection through immersive properties
- Coolcore is the standard cushion cover material that has proven temperature and moisture control. Additional zippered pockets on side of cushion

SIDE POCKETS



1. Special Order Form
 2. CoolCore® is standard on all cushions except the Simplicity G.
 3. Special order sizes available for all cushions except the Glacial SP. Contact sales for additional assistance
 4. Allowances listed are estimations for cushions billed in conjunction with Group 3 power chair bases using the KU modifier and may vary by state