




CB COMBO



HEAD SOLUTIONS COMBO OVERVIEW

The Combo is a triple pad head support that allows movement of the lateral pads from flat (even with the middle pad) to 45°, providing lateral support to the head. The unique contour of the side pads eliminates interference with the ears so that hearing is not impacted and the client's ears do not fold over with head turning. Also, this shape does not interfere with the client's peripheral vision.

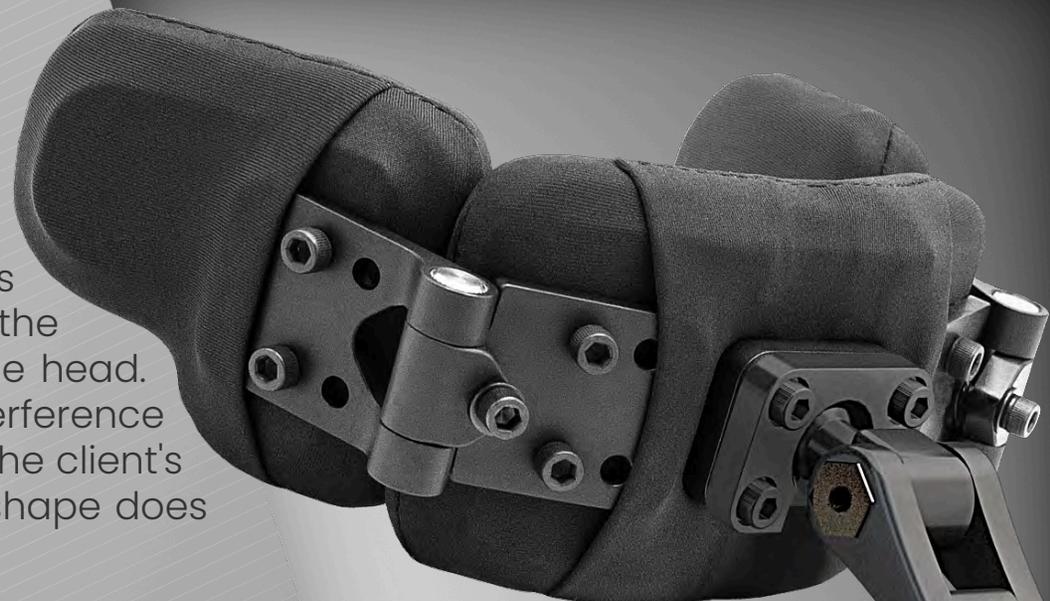

**INDIVIDUALLY
ADJUSTABLE
LATERAL PADS**




**LOW PROFILE
CONFIGURATION**



**STANDARD & CUSTOM
LATERAL PADS**

The standard wedge height can be customized upon request to help prevent ear friction against the lateral wing pads.

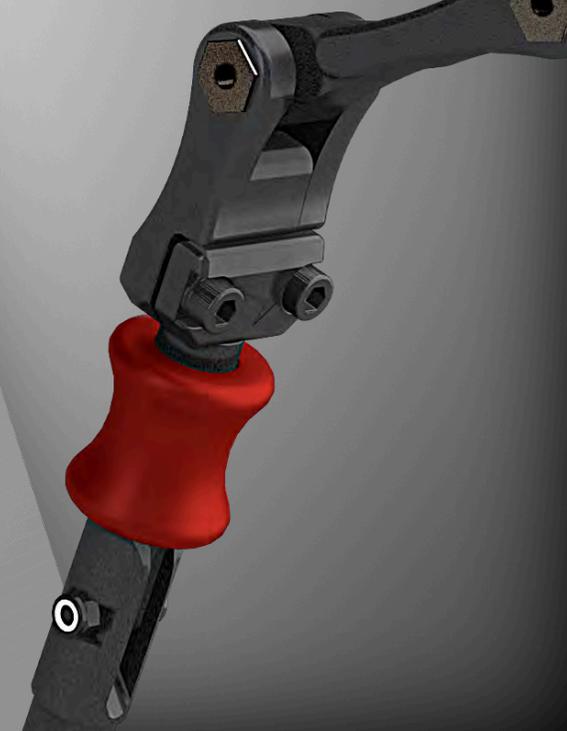

COMBO OPTIONS

SIMPLIFIED OPTIONS




NEED A **LARGER PAD WITH ADJUSTABLE LATERALS?**
CHECK OUT THE **ACP AT THE COMFORT PLUS SERIES**



A



B

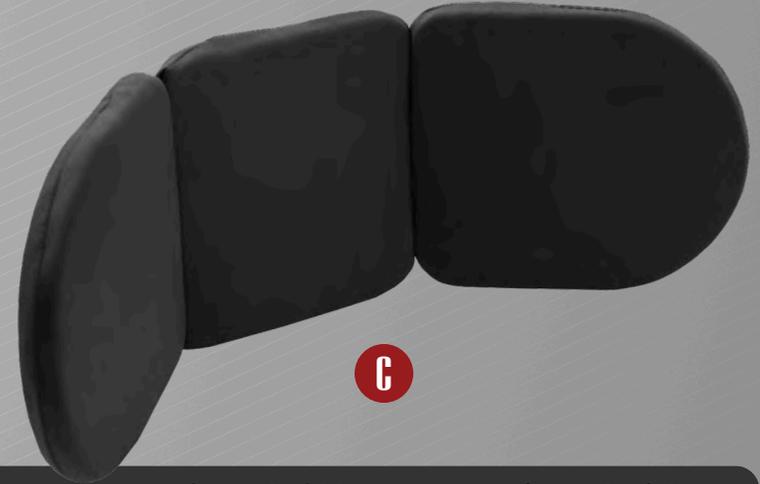

C

	I.N.	Description	Main Pad Width	Main Pad Height	Wing Pad Width	Wing Pad Height
A	CB315M	Metric Small Combo	5in (12.70cm)	3¼in (8.26cm)	3¼in (8.26cm)	3¼in (8.26cm)
	CB315	Small Combo	5in (12.70cm)	3¼in (8.26cm)	3¼in (8.26cm)	3¼in (8.26cm)
B	CB305M	Metric Large Combo	5¾in (14.61cm)	4in (10.16cm)	4¾in (12.07cm)	4in (10.16cm)
	CB305	Large Combo	5¾in (14.61cm)	4in (10.16cm)	4¾in (12.07cm)	4in (10.16cm)
C	CB325M	Metric Extra Large Combo	5¾in (14.61cm)	6¼in (15.88cm)	7in (17.78cm)	6¼in (15.88cm)
	CB325	Extra Large Combo	5¾in (14.61cm)	6¼in (15.88cm)	7in (17.78cm)	6¼in (15.88cm)

Only metric version is compatible with the new UniLink™ head solutions hardware

Adjustable Lateral Hinges



Hinges can extend overall flat pad width **1in (2.54cm)** per side at **1/2in (12.70mm)** increments for a total of **2in (5.08cm)**

Wing Wedge Height



The standard wedge for the small pad is **2in (5.08cm)** and large **2¾in (6.99cm)** but this can be customized at ordering

Removable Pads




5/8in (15.88mm) closed cell latex free foam attached with a **hook & loop** material to the frame for easy removal

COMBO ADD-ONS SWING AWAY

EXTENDS **COMBO**
POSITIONING & **MORE**



MORE **DETAILS** ABOUT
THE **SWING AWAY RECEIVER**
<https://stlpro.site/stealth-swingaway>

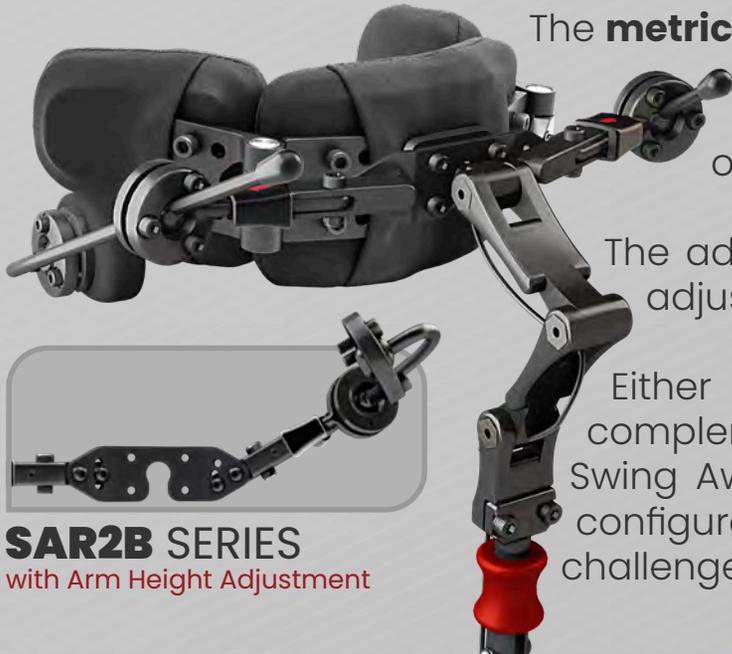

The **metric Combo** is compatible with the **SAR2B** (Low Profile) swing away series while the **original** is compatible with the **SUS9** series. This add-on will allow a caretaker or technician, to position any mini-joystick, switches or temporal pads where they are needed most.



SAR2B SERIES
with Arm Height Adjustment

The advantage of the SAR2B series over the SUS9 series is the arm height adjustment feature which has been introduced with the new lighter design.

Either to enhance positioning or to complement alternative control options, the Swing Away will let a caretaker or technician configure the Combo to offer solutions to any challenge.



SUS9 SERIES

ALTERNATIVE TECHNOLOGY

The **Combo** can be used as a support for tilt/recline. It can supplement with embedded switches on a tray driving control configuration, or any other configuration. It is compatible with our **i-Drive**® Advanced Alternative Control System and **i-Connect**™ sensor switches power router. The **i-Drive Tri-Array** is based on the Combo.

<https://stlpro.site/i-drive>

<https://stlpro.site/stealth-icconnect>

